


Forest Hub

USER GUIDE

Introduction

The Forest Joystick Mouse Hub allows you to connect an analog joystick and up to four (4) assistive switches to use as a USB HID Mouse or USB HID Gamepad. An additional assistive switch can be used to cycle between 3 slots and switch between Mouse mode and Gamepad mode.

Features

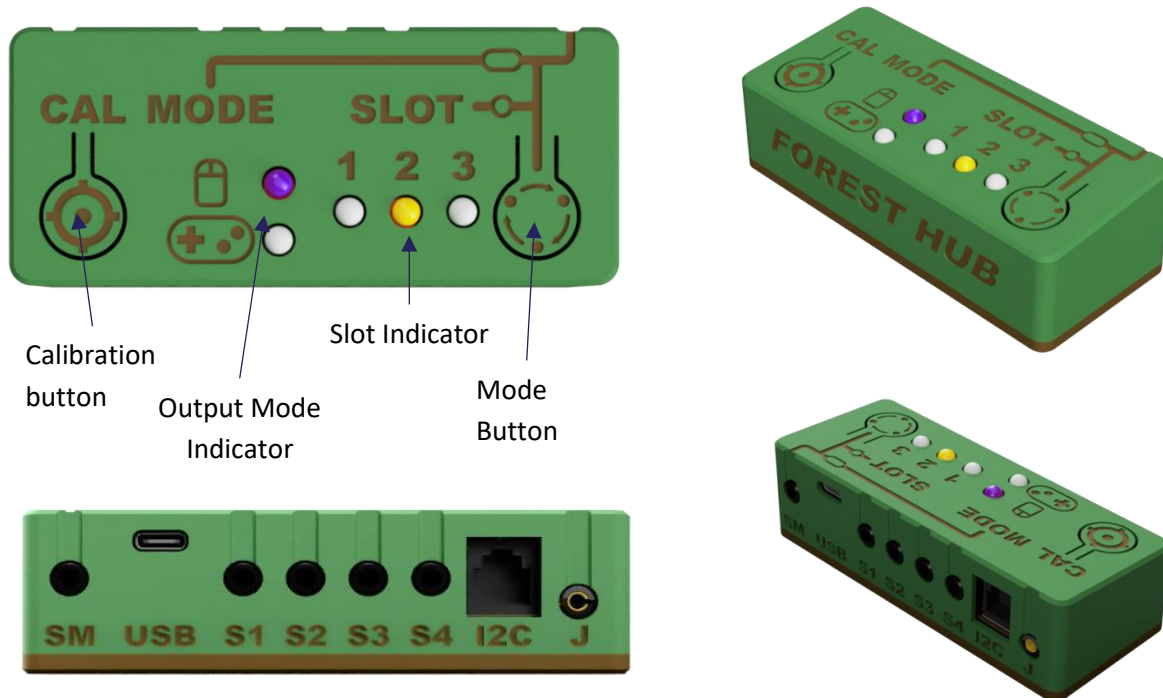

Usage

Setup

1. Connect an analog joystick to the Joystick port (J).
2. Connect assistive switches to the Switch Input ports (S1, S2, S3, S4).
3. Connect an additional assistive switch to the Slot/Mode input port (SM).
4. Connect the USB-C cable to the Hub, then to the Host Device.

Calibration

1. Release the joystick to neutral.
2. Press the calibration button to reset neutral position.

Forest Hub

USER GUIDE

Mouse Mode		Gamepad Mode
<ul style="list-style-type: none"> Move the joystick to control the cursor. Press the S1 button to left click. Press the S2 button to middle click. Press the S3 button to right click. Press the S4 button to toggle scroll. 		<ul style="list-style-type: none"> Move the joystick to control the joystick. Press the S1 button to activate Button 1. Press the S2 button to activate Button 2. Press the S3 button to activate Button 3. Press the S4 button to activate Button 4.
Control	Short Press	Long Press
SW Switch	Cycle Slot	Change mode
Mode Button	Cycle Slot	Change mode
Calibration Button	Neutral position reset	

Host Device Compatibility

- USB HID Mouse: PC, Mac, Android Smartphone, Apple Smartphone (iOS 13+)
- USB HID Gamepad: PC, Xbox Adaptive Controller

Specifications

Enclosure Height	32 mm
Enclosure Depth	46 mm
Enclosure Length	108 mm

Cleaning

The outside of the enclosure can be cleaned with a damp cloth. Do not use hot water. Do not submerge.